



(formerly Warwick Evans Optical Co. Ltd.)

# Keystone VS 6 Medical Vision Screener

*For Medical Practitioners, Transport Companies, Emergency Services etc.*

*Dark Adaptation Exam option keeps your company out of the dark*

## New & Standard VS 6 Medical Features

- Glare Recovery
- Contrast Sensitivity
- Test acuity, phorias, stereopsis
- Enlarged Viewing Lens to Accommodate Multifocals
- Elliptech Soft Touch Membrane Control Panel
- Long-Life True White Light LEDs
- Horizontal and Vertical Peripheral Vision Test
- Night Vision Test
- Visionary™ Software (100 Records)
- 3-Year Warranty
- UK Specification
- No calibration required



### Simple

Simple and easy to use. The VS 6 is layperson friendly requiring minimal examiner training.

### Efficient

The exam process is rapid yet accurate, allowing for maximum efficiency as the screening process averages under 5 minutes per examinee.

### Well-Built

With eight decades of experience, the reliability of Keystone View's engineering and manufacturing ensures your VS 6 Vision Screener will require little or no maintenance.

### Convenient

The VS 6 smartly styled housing provides a small, self-contained, and light-weight vision screener that is easily moved to any location in or out of the office.

### Precise


The VS 6 replicates "real life" vision by using reflective light provided by two long-life true white (5000K) light emitting diodes.

*Eight decades of vision screening leadership*

Sole distributors for UK and Ireland

# VS 6 Medical Vision Screener

The VS 6 Medical blends technological convenience with more than 80 years of vision testing research, expertise and excellence. Like its predecessors, the VS 6 Medical features several unique competitive advantages—such as testing in reflected light rather than with backlit targets to mirror everyday visual function. Now with larger viewing lens the VS 6 accommodates examinees who use progressive lens glasses. The VS 6 is the only screener on the market that comes with glare recovery and contrast sensitivity tests and horizontal and vertical perimeter tests.

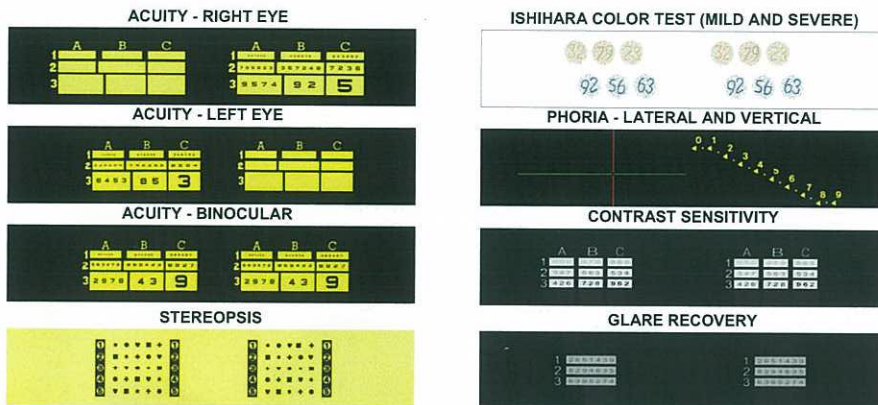

## Available with Visionary Software



- Allows control of the VS 6 test instrument remotely from the computer work station.
- Compatible with most Microsoft Windows based PCs and laptop computers.
- User friendly design makes for “clicking through” the testing procedure.
- Automatically stores results in an easy-to-access client/patient database.
- Multiple reports are easily generated.
- Easily upgrade to unlimited records with additional payment.

*“For recording results, Visionary Testing software lets you customize a database to your specific needs, so you can easily access, analyze and report your results the way you want to see them.”*

Gifford Mast  
President Mast Concepts Inc.



(formerly Warwick Evans Optical Co. Ltd.)

| FEATURE   | BENEFIT   |
|---|---|
| <b>Glare Recovery Test</b>                      | Measures the speed and efficiency by which the examinee re-adapts to night conditions after being flooded with a sudden bright light.   |
| <b>Contrast Sensitivity Exam</b>                | Examinees are presented with nine blocks of numbers varying in contrast from 10 to 90 percent under dusk/night driving conditions. All tests are given at 20/70 (6/21) acuity.  |
| <b>Elliptech Soft-Touch Hand Control</b>        | User Friendly design makes for intuitive control of every test, accessory, and function of your Keystone View VS 6 Vision Screener. Tests can be advanced/reversed automatically.   |
| <b>Unique Long-Life True White LED Lighting</b> | Test targets are illuminated by reflective light the method preferred by most practitioners as it simulates natural visual function and is unaffected by outside factors.   |
| <b>Integrated Head Sensor</b>                   | The integrated head position sensor is utilized to verify correct patient head positioning. This feature is controlled by the Elliptech Soft-Touch Hand Control.  |
| <b>Peripheral Vision Testing</b>                | Standard on all VS 6 Medical vision screeners, peripheral vision testing is performed in the horizontal visual field at 85°, 70°, 55° and nasal 45° for each eye. The vertical field is tested at 35° both above and below the natural line of sight. |

Contact details: NRK Medical Ltd

Email: [sales@nrkmedical.com](mailto:sales@nrkmedical.com) | Web: [www.nrkmedical.com](http://www.nrkmedical.com) | Tel: 020 8888 0051

22 Palace Road, Bounds Green, London N11 2PS

Company number: 430194